


Stryker Prime X X-Ray Stretcher

Image Quality.
Access.
Mobility.

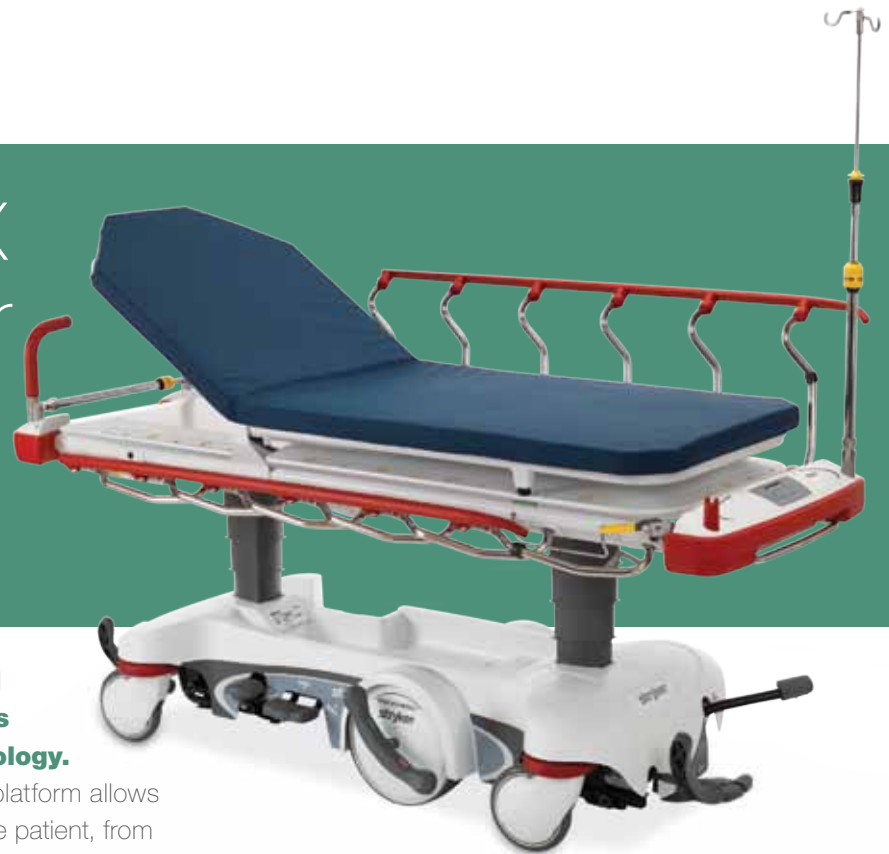

Stryker Prime X offers the safety and efficiency of the Stryker Prime Series stretcher, now with ClearView Technology.

The unique dual-deck design of the patient platform allows positioning of cassettes at any point under the patient, from head to foot and side to side. Open architecture accommodates virtually any size of cassette. This versatile solution helps caregivers deliver exceptional patient care throughout the hospital by combining a mobile platform with superior image quality.



ClearView Technology combined with a dual-deck design offers the image quality you would expect from an x-ray table with the convenience of a mobile platform.



Prime X offers advanced mobility solutions such as Big Wheel, which reduces push force and enhances manoeuvrability, or the Zoom Motorized Drive System,* which virtually eliminates manual pushing.



Prime X offers 360 degrees of access to quickly and easily place the cassette for full body imaging.

Dual-deck design allows cassette to be adjusted to the patient.

Open architecture accommodates virtually any size of cassette.

Bevelled edges and visual alignment guides on the cassette platform help make cassette placement quick and easy.

Electric lift makes raising the stretcher easy by simply stepping on the pedal.

700 lb./317 kg. weight capacity accommodates your heavier patients.

Articulating backrest and foot section allows complete access to the dual deck for easy cleaning.

Glideaway Siderails provide exceptional coverage and tuck neatly under the litter for unobstructed lateral imaging and storage when you don't need them.

*Pending 510(k) approval, not yet available for sale

Stryker Prime X X-Ray Stretcher

Stryker SA
Cité-Centre
Grand-Rue 90
1820 Montreux
Switzerland

t : +41 21 966 12 01
f : +41 21 966 12 00

www.stryker.eu



1105X

Standard Features

- ClearView Technology
- Thermoformed ABS base with extra capacity storage tray
- 700 lb. weight capacity
- 30-inch wide patient surface with 3-inch Enhanced Comfort Mattress
- Litter-mounted oxygen bottle holder
- Integrated transfer board
- Power-washable design (1105X and 1115X only)
- Four-wheel, center-locking brake system with dual end activators
- Dual, side-mounted foot control hydraulics with uni-lower pedal
- Glideaway Siderails with integrated foot end push/pull handles
- 4 integrated IV receptacles
- +/- Trendelenburg
- 8-inch Omni Surface Casters with integrated wheel covers
- Pneumatic backrest/stationary foot
- Constant descent hydraulics
- Dual end brake/steer controls
- 2-year warranty

1105X: All Standard Features Plus

- Fifth Wheel steering

1115X: All Standard Features Plus

- Big Wheel steering

1125X: All Standard Features Plus

- Zoom Motorized Drive System
[Pending 510(k) approval, not yet available for sale]

Optional Features

- 4-inch (10 cm) Enhanced Comfort Mattress
- 4-inch (10 cm) Ultra Comfort Mattress
- Four-sided brake and steer controls (1105X and 1115X only)
- Pop-up push handles (1105X and 1115X only)
- Three-sided hydraulic controls (1105X and 1115X only)
- Defibrillator tray
- Defibrillator tray/foot extender/charting surface
- Dual-ended siderail release
- Footboard/chartholder
- IV caddy
- IV poles (permanent, folding two- or three-stage)
- Serving tray
- Serving tray holder/footboard
- Siderail pads
- Upright oxygen bottle holder
- Restraint straps
- Powered lift (1105X and 1115X only)

This document is intended solely for the use of healthcare professionals.

A surgeon must always rely on his or her own professional clinical judgment when deciding whether to use a particular product when treating a particular patient. Stryker does not dispense medical advice and recommends that surgeons be trained in the use of any particular product before using it in surgery.

The information presented is intended to demonstrate the breadth of Stryker product offerings. A surgeon must always refer to the package insert, product label and/or instructions for use before using any Stryker product.

Products may not be available in all markets because product availability is subject to the regulatory and/or medical practices in individual markets. Please contact your Stryker representative if you have questions about the availability of Stryker products in your area.

Stryker Corporation or its divisions or other corporate affiliated entities own, use or have applied for the following trademarks or service marks: **Stryker, Big Wheel, ClearView, Prime X & Zoom**. All other trademarks are trademarks of their respective owners or holders.

The products listed above are CE marked according to the Medical Device Directive 93/42/EEC.



1115X

Specifications

Model Number	1105X (Prime X with Fifth Wheel)* 1115X (Prime X with Big Wheel)* 1125X (Prime X with Zoom Drive)*
Overall Length	
(1105X and 1115X)	85" (216 cm)
(1125X)	86" (218 cm)
Overall Width	
Siderails Up	38" (96.5 cm)
Siderails Down (1105X and 1115X)	30.5" (77.5 cm)
Siderails Down (1125X)	31" (78.7 cm)
Patient Surface	30" x 75.25" (76 x 191 cm)
Weight Capacity	700 lb. (317 kg)
Height Range	
High	37.25" (94.6 cm)
Low (1105X and 1115X)	24" (61 cm)
Low (1125X)	26.25" (66.68 cm)
Litter Positioning	
Backrest	0°– 90°
Foot Section	0°– 40°
Trend./Reverse Trend. (1105X and 1115X)	±17°
Trend./Reverse Trend. (1125X)	±16°
Siderails	14" x 57.75" (36 x 147 cm)
Caster Diameter	8" (20 cm)
EMC	IEC 60601-1-2

Warranty

Two years parts, labor and travel

Stryker reserves the right to change specifications without notice.



Mkt Lit-MTX626121311EN Rev A
MTX-RRD/GS 02/12
Copyright © 2012 Stryker

